

mooC



# CONTENTS

1. Hydrafrac
2. Picode
3. Picode Lite
4. C-PLus
5. Code Factor
6. Em-Code
7. CoolShape
8. Em-Cool
9. Excimer
10. Nodd Pro-V
11. Meso Skinjector
12. Facecode
13. Visage
14. Wavelength Pro-X

# ABOUT CODE AESTHETICS

## INTRODUCTION

---

Code aesthetics is a professional distributor of aesthetics & medical laser equipments in india, established in July 2008. Code aesthetics has its own marketing & development center, studios center, sales and after-sales departments, and offer professional technology, support and clinic data. Code aesthetics team can communicate fluently in Hindi, English, with the customers.

Our Mission: Create a world-renowned brand, build a first-class team, and win the glory for india

CODE Aesthetics offers an extensive array of laser machines, each designed for a specific treatment and indication, together covering a full spectrum of solutions.



# hydrfrac

Redefined Hydration

HydraFrac is a dermal micro-needling machine that offers simultaneous delivery of infused ingredients. It makes the treatment a quick, pain-free and effective treatment with no downtime.

It is the latest innovative platform merged with all the technologies.

INTRODUCING  
NEW ERA OF MEDI-FACIALS  
HYDRAFRAC...

- Acne scarring
- Skin dehydration
- Superficial wrinkles & Fine lines
- Uneven skin tone, & Pigmentation

HydraFrac is (USFDA) approved platform. It involves Radio Frequency (RF), Hydra Derma Abrasion (HDA), Vacuum Assisted derma roller (VADR), SONO, electro-ion and Glow-O to achieve remarkable results.

HydraFrac is the best option for skin Rejuvenating, Hydration, Infusion, and pigmentation.

The HydraFrac machine takes micro-needling to a whole new level.

## TECHNICAL SPECIFICATION

✔ Ultrasound  
60 KHz/2ww/cm

✔ Screen  
10inch tft

✔ Input Voltage  
110V/220V 50/60Hz

✔ Radio Frequency  
2 MHz, Multi Polar (6pol)

✔ Vacuum Range  
Max 680nm HG

✔ Noise Level  
45DB

✔ Ion lifting  
500Hz (Digital Ion Lifting)



Electro-Ion



HydraFrac Sono



HydraFrac MRF



GLOW-O



Aquapeel



### Exfoliator

Gently Exfoliates Dead skin cells

Active ingredients  
Glycolic acid



### Controller

Cleanse, oil control & open pores

Active Ingredients  
Salicylic acid



### Accelerator

Accelerate the cell metabolism

Active ingredients  
EGF & galactomyces



### Cleaner

Aquapeel handpiece cleaner

Skin Rejuvenation Gel  
6x0.68fl.oz.20ml.e

Mild Exfoliation Capsule  
6x0.01 lb 6gr.e




## HYDRAROLLER

Non-Invasive Microneedling Apparatus



PATENTED DESIGN

Hydraroller is a non-invasive micro protrusion attachment consisting of our unique tightly placed micro protrusions which debris the skin surface safely and controllably.

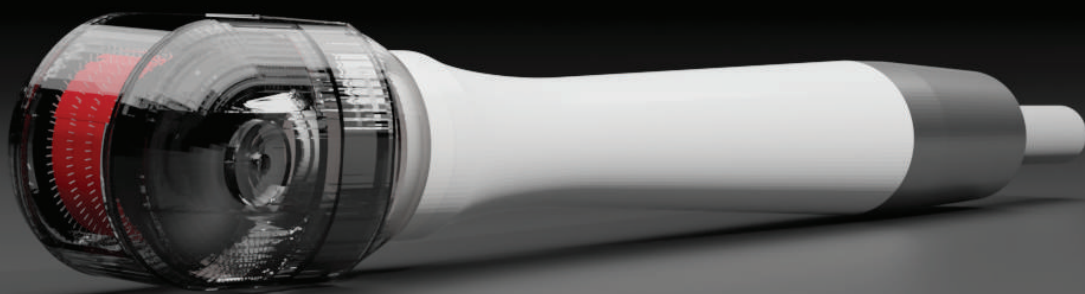

### Hydra Needle Tip

Titanium needles best suited for medical micro needling treatment at DE junction (0.25mm)

- More effective NON-invasive infusion - for better results and serum permeation with each pass.
- Vacuum Assisted-more of the serum is absorbed.

# PI-CODE

PICO SECOND LASER



Chloasma  
Removal



Tattoo  
Removal



Ota Nevus  
Bluemole etc.



Remove  
Spots

## FEATURES

- Shortest pulse duration, enabling high peak power
- Unique double pulse mode
- Supreme feature set
- Improved ROI with excellent results in fewer treatments
- A versatile platform enabling multiple patient-tailored treatment options
- Suitable for all skin types

## TECHNICAL SPECIFICATION

LASER TYPE	PI-Code ND-YAG Laser
WAVELENGTH	1064nm, 532nm
FREQUENCY	1 - 10Hz
PULSE WIDTH	750PS
LASER OUTPUT MODE	Flat top
ENERGY	0.3 - 31.8J/cm <sup>2</sup>
SPOT SIZE	2 - 10mm
COOLING SYSTEM	Air cooling + water cooling
VOLTAGE	220V/110V
DIMENSION	112cm*47.5cm*102cm
WEIGHT	N.W. 109.6KG G.W. 164.6KG

AD55

## INTRODUCTION

PI-Code presents a new era in skin rejuvenation and tattoo removal, offering optimal results for the wrinkle reduction, removal of colorful tattoos, and benign pigmented lesions with just a few treatments. PI-Code mode of operation is based on delivering ultra-short picoseconds pulses of energy into the skin tissue.

## ACCESSORIES

PI-CODE offers an extensive array of applicators, each designed for a specific treatment and indication, together covering a full spectrum of solutions.



### 1. PI-Code Handpiece

Offers 9 spot sizes options, ranging between 2 and 10mm, addressing various types, severities, and depths of pigmented lesions and tattoos. The applicator can be used both in Nd: YAG (1064nm) mode and Double Nd: YAG (532nm) mode, Using a simple on-screen setup.



### 2. PI-CODE HANDPIECE

Employs a fractional delivery method that creates pixel-sized perforations in a 7x7 Non-ablative pattern, leaving the surrounding tissue intact. These micro injury sites trigger a wound healing process that strengthens collagen and stimulates new collagen, completely rejuvenating the target tissue. Five distinct treatment depths are available for maximum flexibility and precision.



# PI-CODE LITE

PICO SECOND LASER

## INTRODUCTION

Pi-Code presents a new era in skin rejuvenation and tattoo removal, offering optimal results for the wrinkle reduction, removal of colorful tattoos, and benign pigmented lesions with just a few treatments.

Pi-Code mode of operation is based on delivering ultra-short picoseconds pulses of energy into the skin tissue.

### SKIN REJUVENATION

Skin rejuvenation treatment protocols are based on fractional energy delivery at various skin tissue depths.

### TATTOO REMOVAL

Pi-Code mode of operation is based on delivering ultra-short picoseconds pulses of energy into the skin tissue.

### PIGMENTED LESIONS

Pi-Code is able to treat both superficial & deep pigmented lesions. The dual-wavelength option allows us to treat epidermal-dermal lesions with its 532nm wavelength and deep pigmented lesions with The 1064nm to target melanin.



## TECHNICAL SPECIFICATION

LASER TYPE	PI-Code ND-YAG Laser
WAVELENGTH	1064nm, 532nm
FREQUENCY	1 - 10Hz
PULSE WIDTH	750PS
LASER OUTPUT MODE	Flat top
ENERGY	0.3 - 31.8J/cm <sup>2</sup>
SPOT SIZE	2 - 10mm
COOLING SYSTEM	Air cooling + water cooling
VOLTAGE	220V/110V
DIMENSION	112cm*47.5cm*102cm
WEIGHT	N.W. 109.6KG G.W. 164.6KG



## ACCESSORIES

PI-CODE offers an extensive array of applicators, each designed for a specific treatment and indication, together covering a full spectrum of solutions.

### 1. PI-CODE HANDPIECE

Offers 9 spot sizes options, ranging between 2 and 10mm, addressing various types, severities, and depths of pigmented lesions and tattoos. The applicator can be used both in Nd: YAG (1064nm) mode and Double Nd: YAG (532nm) mode, Using a simple on-screen setup.

### 2. PI-CODE HANDPIECE

Employs a fractional delivery method that creates pixel-sized perforations in a 7×7 Non-ablative pattern, leaving the surrounding tissue intact. These micro injury sites trigger a wound healing process that strengthens collagen and stimulates new collagen, completely rejuvenating the target tissue. Five distinct treatment depths are available for maximum flexibility and precision.

## THERAPY PRINCIPLE

### Therapy of epidermis and dermis pigmentation

Utilizing the explosive effect of the ND: YAG laser, the laser permeate the epidermis into the dermis which includes amount of pigment mass. Since the laser pulses in nanosecond but with super high energy, the shot pigments mass swell quickly and breaks into small pieces, which will be eliminated through metabolic system.

### Therapy of the dilatation of capillary vessel

Utilizing the heat effect of the ND:Yag Laser, the laser is absorbed by the monochrome in the capillary vessel, Then the capillary vessel become blocked the disappear.

## TREATMENTS

Wrinkles Reduction

Pigmented Lesions

Colored Ink Tattoos




Skin Rejuvenation



Pigmentation



Dark Ink Tattoos

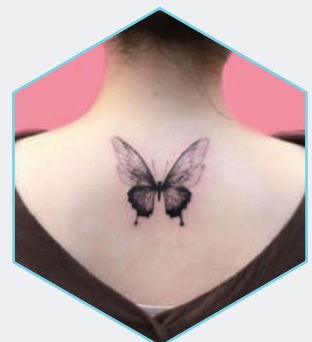

## Technical specifications

### Basic operation

Type of laser equipment	Q - Switch ND:YAG Laser
Lens type	1064-F50, 1064F-FO, 532F-5O

Power supply	220-240V, 50-60Hz
Frequency	1 - 10Hz
Energy	1000 mj

### Technical characteristics

Mode	1064nm, 532nm
Pulse width	10ns
Cooling system	Water cooling , Air cooling


Weight	25kg
--------	------

CODE ESTHETICS

# C-PLUS

Q-SWITCHED ND:YAG LASER

C plus+ presents the powerful triple mode nd: yag laser treatment solution available, featuring q-switched, laser toning, and carbon laser on a single, dedicated platform. With double pulse technology offering up to 1000mj per pulse, an extraordinarily short pulse width of 7~12 nanoseconds, unique technological innovations including depth control and fractionated delivery. C- plus+ is unmatched in its power and versatility for treating a wide range of aesthetic dermal indications.

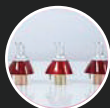

## WHAT IS C-PLUS ?

Principle of q-switched nd:yag laser is that it generates strong energy in a very short time to instantly destroy pigments lesions in both epidermis and dermis without thermal damage to the epidermis and surrounding skin, follows by macrophages absorb destroyed pigment lesion in the human body



Selective Reaction  
Only For Pigment




Pigment Removal



## ACCESSORIES

### MAIN APPLICATOR

Treats pigmented lesions and tattoos using a homogenized laser beam profile with uniform energy intensity, preventing hot spots. A square beam with  $3 \times 3 \text{ mm}^2$  or  $5 \times 5 \text{ mm}^2$  spot sizes allow for coverage of treatment areas without overlap. The homogenizer handpiece may be used in q-switched laser mode

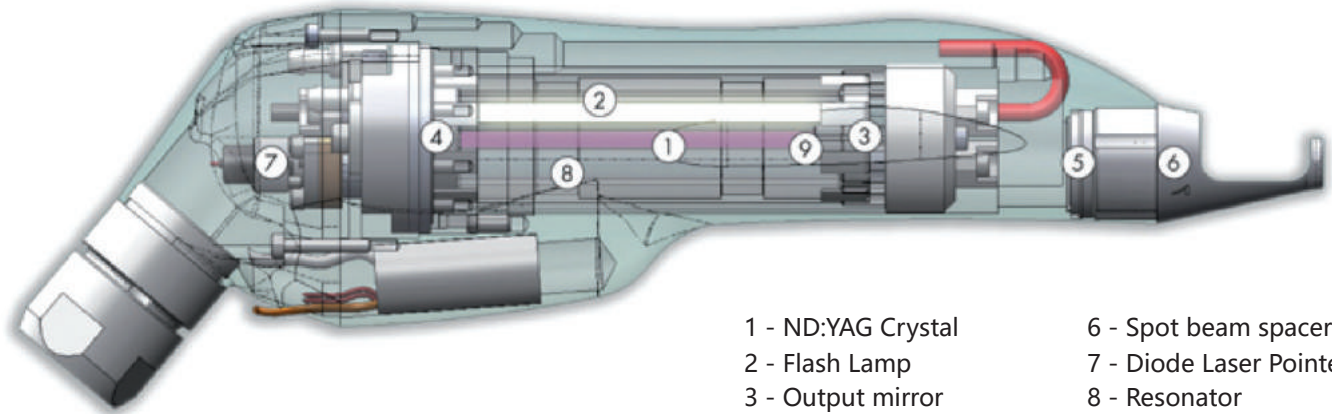

- |                         |                         |
|-------------------------|-------------------------|
| 1 - ND:YAG Crystal      | 6 - Spot beam spacer    |
| 2 - Flash Lamp          | 7 - Diode Laser Pointer |
| 3 - Output mirror       | 8 - Resonator           |
| 4 - Reflective mirror   | 9 - Q-Switched Crystal  |
| 5 - Beam Expansion lens |                         |

### APPLICATOR HEADS

The 3 handpiece delivers parallel beams of energy to the target tissue with minimal dispersion, regardless of distance from the skin. This allows practitioners to administer treatment without having to maintain continuous contact with the skin while also allowing better visualization of larger treatment areas. With a 5mm spot size, this handpiece offers an excellent coverage rate and high-speed treatment.



1064NM-F0




#### USING 1064NM- F0

After applying the carbon cream the carbon particles absorb the laser light and explode instantly, removing impurities and keratin on the skin surface and helping to regenerate collagen



1064NM-F50

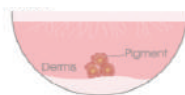


#### USING 1064NM- F50

For pigment disease at dermis layer with strong energy tattoo, Nevus of OTA



532NM-F50

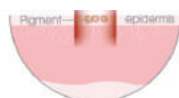

#### USING 532NM -F50

For pigment disease at epidermis layer with strong energy freckles, age spot.

# CO<sub>2</sub> FRACTIONAL LASER

## INTRODUCTION

Factor code is a state-of-the-art, minimally-invasive outpatient solution that incorporates CO<sub>2</sub> technology to provide optimal care for various feminine concerns. The procedure is fast, painless, and completely safe, and provides highly effective results. It's a versatile solution that enables medical practitioners to treat various feminine indications quickly and effectively. Factor code opens up your clinic to a wide range of patients, varying in age and stage in life, broadening your reach – and adding real value to your practice.



Fractional CO<sub>2</sub> Laser For Vaginal Tightening



Fractional CO<sub>2</sub> Laser For Vulva Treatment



CO<sub>2</sub> For Surgery



Fractional CO<sub>2</sub> Laser For Facial Treatment

CODE  
FACTOR

## TECHNICAL SPECIFICATION

WAVELENGTH	10.6um
WORKING MODES	Continuous/Ultra-pulse/non-sequential Pulse
POWER	40W
PULSE WIDTH	0.1-10ms adjustable Step 0.1ms
LASER	DC Closed RF Laser Generator
INDICATOR LIGHT	635nm infrared
LASER ORIGIN	USA
SCANNING GRAPHICS	Square, Rectangle, Circle, Diamond, Triangle, Ellipse.
LANGUAGE	English
LCD DISPLAY	15 Inch 16- Bit color LCD Display
COOLING SYSTEM	AIR
VOLTAGE	110v-240 AC
SIZE	530*570*1380mm
WEIGHT	55kg

HAND PIECE IDENTIFICATION AUTOMATICALLY  
 ERGONOMIC DESIGN, EASY TO OPERATE



**GYNECOLOGY HANPIECES**  
 ( F127MM )

Rebuiled vaginal mucosa,  
 Tighten vagina and improve  
 The vulva color

01.

**Gynecology Handpieces (f127mm)**

Rebuild vaginal mucosa,  
 Tighten vagina and improve the vulva color.

The Vaginal handpiece is a unique, single-use hygienic probe, uniquely designed to fit varying vaginal sizes and ensure maximum patient comfort. During treatment, the physician rotates the probe in a 360-degree pattern to cover the entire vaginal wall, providing optimal coverage and efficiency.



**LIGHT GUIDE ARM**

Golden standard of anti aging treatment  
 acne and scar improvement

02.

**Scan Handpieces**

Versatile fractional scanner for efficient ablative procedures,  
 ideal for skin texture and pigmentation

LiteScan is a handheld microprocessor-controlled laser scanner for char-free ablative procedures that allows for rapid coverage of large areas. It is designed for use in a variety of dermatological, aesthetic and surgical applications. Multiple scanning options are available including square or circle, in a spiral, grid pattern as well as straight or curved lines.



**SURGICAL HANDPIECE**  
 ( F50MM F100MM OPTIONAL )

Accurate energy  
 Skin tissue instant gasification

03.

**Surgical Handpieces..**

(f50mm,f100mm optional)

This unique technology approach enables the use of much higher energy levels via longer pulse durations, creating a greater coagulation and thermal effect with minimal ablation and significantly reduced patient downtime.

The deep thermal effect of the Pixel applicator makes it ideal for skin rejuvenation, treating fine lines and wrinkles.



**MODES OF OPERATION...**

**THE CO<sub>2</sub> LASER SYSTEM OFFERS 4 MODES OF ENERGY DELIVERY**

1. CW MODE (CONTINUOUS WAVE): \_\_\_\_\_  
 Emits one continuous and consistent laser energy, typically used for excision or cutting.

2. PULSE MODE (SINGLE-PULSE): \_\_\_\_\_  
 Produces a single pulse of laser energy for more supervised and accurate energy delivery.

3. REPEAT MODE: \_\_\_\_\_  
 Produces a series of either short or long pulses of laser energy at variable speeds.

4. SUPER PULSE MODE: \_\_\_\_\_  
 Super pulse mode provides the highest peak power with the shortest pulse duration, Enabling the deepest penetration with minimal residual thermal damage. This mode is designed for performing ablation without inflammation and achieves highly effective results with minimal downtime and enhanced patient comfort.

## BODY CONTOURING



# EM-CODE

BUILD MUSCLE & BURN FAT



NON-INVASIVE



BUTTOCK PROCEDURE



NO DOWNTIME



NO PAIN



NON RADIATION




BODY SCULPTING

## INTRODUCTION

EM-Code is an effective way to improve your fitness level, body image and well-being in far less time.

With latest generation technology, high-tech non-water cooling system, high quality insulated material, scientific magnetic pulses, automatic control system, one size is suitable for multiple areas. EM-Code body workout activates over 90% of the muscles fibers at the same time and also deeper muscles are trained

- User Friendly Screen
  - High resolution interface
  - Free sliding operation
  - 45° Adjustable for best viewing
  - Auto adjustment energy range when changing frequency
  - Smart control of total tesla per session
  - Safety low frequency waves without adverse effect

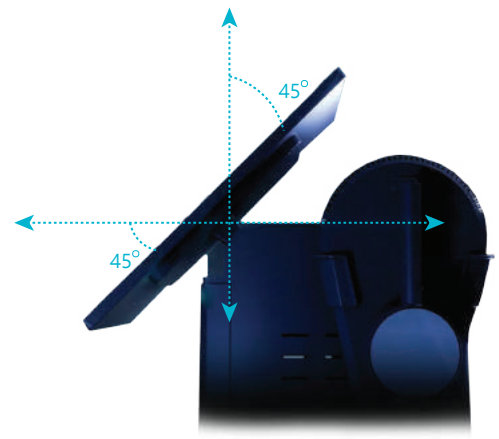

### 4 Working Modes


Intelligent operation, time adjustable from 0 to 60 minutes; energy adjustable from 0 to 100%

- Simple 1 - For people who seldom exercise with less muscle and more fatness  
Default frequencies are 30Hz 38Hz 5Hz
- Simple 2 - For people who do regular exercise with more muscle and less fatness.  
Default frequencies are 20Hz 29Hz 4Hz
- Specialty 1- For professional practitioners to customize protocol A  
The frequencies of F1, F2, F3 are adjustable, and the energy adjustment range is automatically optimized according to the frequency
- Specialty 2 - For professional practitioners to customize protocol B  
The frequencies of F1, F2, F3 are adjustable, and the energy adjustment range is automatically optimized according to the frequency



## Technical Specification

Pulse frequency	1-150Hz
EM Energy	0-100%
RF Energy	0-100%
RF Frequency	6MHZ
Pulse Time	285us
Magnetic Field strength	0.4 - 2.5T
Treatment Time	Up to 60 min
Equipment weight	40kg
User Interface	Touch screen
Operation type	Continuous

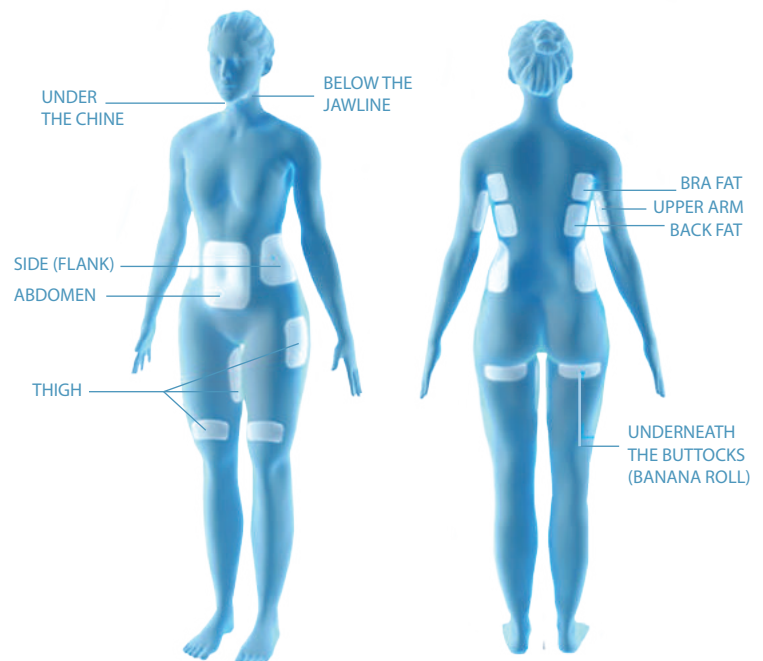

# CoolSHAPE

7<sup>th</sup> GEN CRYOLIPOLYSIS

The latest in non-invasive full body contouring.

## INTRODUCTION

Coolshape uses state -of-the-art freezing technology removes excess subcutaneous fat from a specific part of the body by dissolving fat cells using cool temperature without damaging the surrounding tissue. Once cold, fat cells will decompose naturally and gently be eliminated through the body's normal metabolic processes. Designed with 360° surrounding cooling technology, the multiple cooling cups fit all curves and contours from as high as your chin all the way down to the knees to deliver absolute, non-invasive body contouring your way.



### 7 IN 1 MULTI-SIZE APPLICATOR

Ergonomic design, skin-friendly fit the body arc meet the whole body fat reduction needs from the face to the legs

### Medical flexible LSR silicone cups

Whole shaped molding 360 cups dedicated and soft, comprehensive care for the skin.

### 360 Cooling

Better Cooling effect and less treatment time.

### Liquid Silicone

Soft and tender touch like a baby skin.

### Dual Temperature

Exclusive patent assures safety and efficacy.




## AIM OF THE TREATMENT

Coolshape treatment eliminates fat tissue in numerous areas allowing full body contouring. Mainly, it is a pain-free, non-invasive and definitive treatment for the elimination of fat cells.




## TREATMENT APPLICATION


Unique Patented Technology Ensure Best Results, comfortability and safety

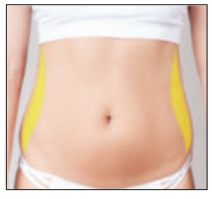


You may struggle with stubborn fat even after jogging and dieting.

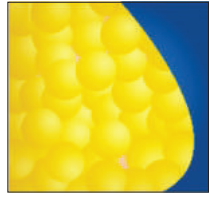


The coolslimming suction your skin and subcutaneous layer.

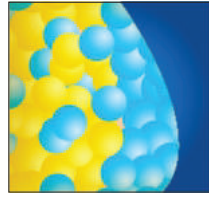


Cooling energy then freezes fat deposit to induce apoptosis.

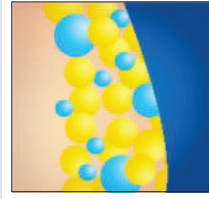

Crystallized fat cells die off through the body's natural metabolism.

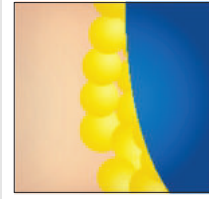

Remaining fat in the e target area became leaner and less dense.

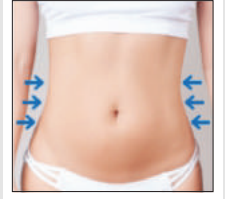

Combine with your healthy regimen to look slimmer and confident.

## Applicators working in unison


Multiple combinations in just one session.

ALL AREAS

ALL ANGLES

ALL TYPES OF TISSUE

ALL PATIENTS

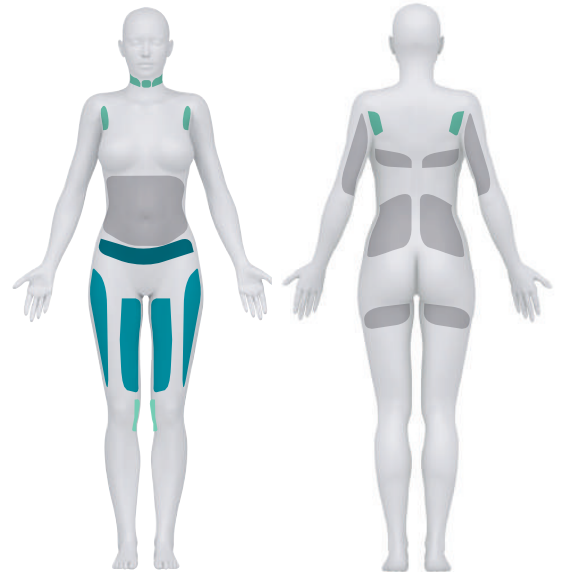

## EFFECTIVE FAT REDUCTION

The newly added RF energy promotes blood circulation, heats the cells in the deep layers of the skin, reduces the thickness of the subcutaneous fat layer, and makes the skin contour plump.





# EM-COOL<sup>HYBRID</sup> C O N T O U R I N G

## INTRODUCTION

**EM-COOL Hybrid** uses state -of-the-art freezing technology removes excesss subcutaneous fat from a specific part of the body by dissolving fat cells using cool temprature without damaging the surrounding tissue. Once cold, fat cells will decompose natural and gently be eliminated through the body's normal metabolic processes. designed with 360. surrounding cooling technology. the multiple cooling cups fit all curves and contours from as high as your chin all the wqay down tothe knees to deliver absolute, non-invasive body contouring your way.

## APPLICATOR INTRODUCTION

The Four- Generation applicators are upgraded with exclusive patent, creating a core competitive advantage the is unique to **EM-Cool hybrid**

## FAST AND EASY APPLICATOR EXCHANGE

With a simple press-and-release mechanism, the EM-Cool Hybrid allows maximum Convenience and makes switching cooling cups as simple as possible. Detaching Cables or turning the system off is not requires even during treatments.

The part description is for reference only, please select the corresponding handle according to the individual body for actual operation.



**Wing Mini**

Double Chin,  
Axilla, Knees



**Flat Mini**

Double Chin,  
Axilla, Knees



**Flat**

Abdomen,Arms,  
Inner Thighs



**Arc-T**

Outer Thighs,  
Inner Thighs



**Wing**

Abdomen,Flanks,  
Love Handles



**Wing Plus**

Abdomen,Bra Line,  
Flanks,Love Handles,  
Below Buttocks




**Wing Max**

Abdomen,Flanks,  
Love Handles

## EFFECTIVE FAT REDUCTION

The newly added RF energy promotes blood circulation, heats the cells in the deep layers of the skin, reduces the thickness of the subcutaneous fat layer, and makes the skin contour plump.

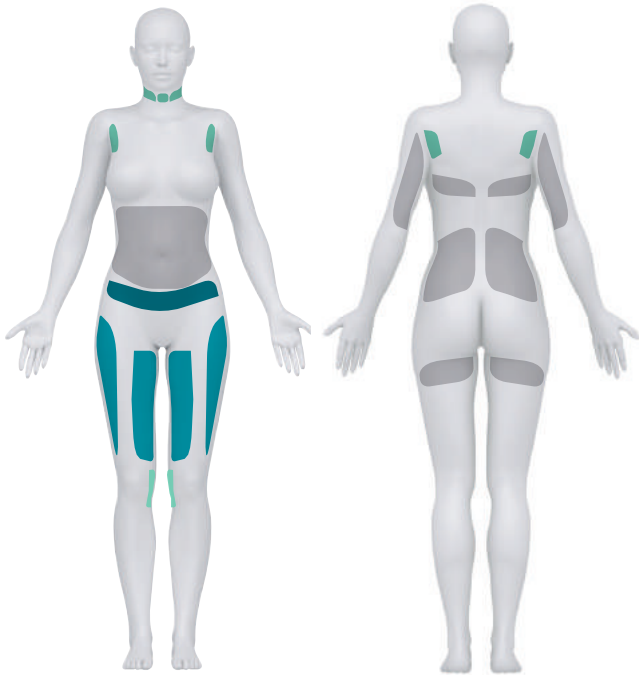

## Applicators working in unison

Multiple combinations  
in just one session.

ALL AREAS

ALL ANGLES

ALL TYPES OF TISSUE

ALL PATIENTS



NON-INVASIVE



BUTTOCK PROCEDURE



NO DOWNTIME



NO PAIN



NON RADIATION



BODY SCULPTING

## \* Technical specifications

### Basic operation

Type of equipment Controlled cooling device for medical use

### Technical characteristics

Minimum temperature - 11°C  
hand-pieces Straight HP, Curved HP, Tight HP  
Double HP, Tiny HP, Oval HP

Power supply 100-240 V  
Maximum power 1000 W  
Cooling capacity in each handpiece 90 W programming 0 °C  
with a Tamb of 25 °C

### Regulation


Class of equipment as standard Class I, type BF  
classification according to directive Class II a  
Marking CE0051  
The device complies with the requirements Directive 2007/47 CE (93/42/CE)  
Of directives (Medical Device)  
Directive 2002/95 CEE (RoHS)  
Directive 2002/96 ECC (RAEE)

dimensions 540 x 440 x 550 mm

weight 39 kg

screen touch screen 10,4"

CERTIFIED BY

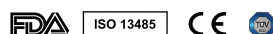

# EXCIMER

## CODE 308 SYSTEM

UVB THERAPY & MONOCHROMATIC THERAPY WITH 308nm



### Code 308 System

#### Monochromatic UVB light therapy for Psoriasis and Vitiligo skin diseases.

The code 308 system is a handheld system for intense monochromatic uvb therapy, treating autoimmune diseases of the skin such as psoriasis and vitiligo. Using a 308nm wavelength which replicates the wavelength of the excimer laser, code 308 system allows you to treat patients quickly and effectively with the benefits of a lighter weight system that offers greater flexibility, lower investment costs, no consumables and greater ease of use.



- Clinically-proven with effectiveness matching that of an excimer laser
- Focused treatment with maximum protection of surrounding skin
- Pain-free with no downtime
- Compact, portable system
- Cost-effective with low operation costs and maintenance-free technology

### EXCIMER SYSTEM VS. EXCIMER LASER

DEVICE	EXCIMER SYSTEM	EXCIMER LASER
Wavelength	308nm	308nm
Clinical Results	Indicated	Identical
Acquisition Cost	Low	High
Maintenance Cost	Minimum	High
Weight	Mobile and handheld system	60kg



## TREATMENTS

The 308nm wavelength is ideally suited for treatment of UV sensitive skin diseases such as psoriasis and vitiligo. The code 308 system may also be used to treat other skin diseases that stem from immune system dysfunction. UVB radiation has an immune modeling effect which can normalize pathologically altered cellular activities. The ideal indications are persistent, localized lesions.

## INDICATIONS



EXCIMER LASER TREATMENT CAN BE USED IN

- Vitiligo: localised, few lesions.
- Psoriasis: localised, palmoplantar, scalp psoriasis.
- Alopecia areata
- Cutaneous T cell lymphoma
- Localised scleroderma
- Granuma annulare

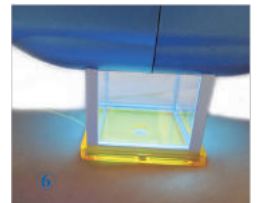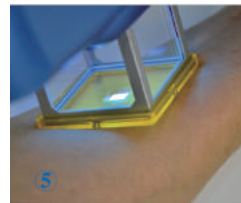

## TECHNOLOGY

The code 308 system uses a 308nm wavelength to effectively target skin lesions. The 16 cm<sup>2</sup> spot can be applied directly into localized plaques, or masks may be mounted to adapt the spot size and shape to smaller areas. A variable power density ranges from 50 to 6,000 mj/cm<sup>2</sup> and can be adapted to The needs of the individual patient. The system enables you to apply high, therapeutically-effective doses of monochromatic UVB light while leaving healthy skin unexposed and protected.

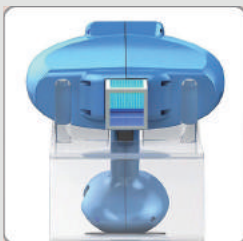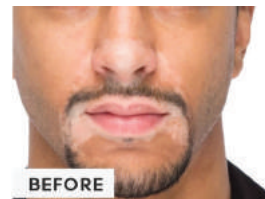

## PRODUCT FEATURE

1. The 308nm Single Wavelength excimer light output via the medium XECI
2. High dose, highly effective and fast
3. Targeted treatment, Maximum protection of healthy skin
4. Up to 16cm<sup>2</sup> spot sized output windows to meet different body parts treatment demand in clinical practice
5. Multiple spot size for short treatment times of larger areas, Fast Treatment
6. No Replacement consumables, Low Maintenance cost
7. User-friendly interface, Convenient and easy to operate
8. Handheld Designed excimer system mobile for therapy



## FRACTIONAL RF APPLICATOR

FRACTIONAL RF applicator for non-invasive wrinkle reduction, suit the patients fear of needle, no down time, recovery soon.

TECHNICAL SPECIFICATION	
System Type	Mono & Bipolar RF
	Bipolar RF
Frequency	2MHz
Total RF treatment Level	Upto Level 10
Dual Operation Mode	Fractional RF Micro-needle
	Fractional RF Sub-Lative
Max.Engery	50W
Disposable Tip	Micro-electrode needle
	Micro-electrode tip
	0.5 - 3 mm (Adjustable)
Display	10.4 inch LCD
Switch	Foot Switch
Dimension(mm)	280mm(W) x 505mm(L) x 1100mm(H)
Weight	20Kg
Electrical Requirements	AC/100-240 V, 50-60 Hz

## INTRODUCTION

NODD Fractional RF Microneedle is designed specially for delivery RF energy directly to epidermis or dermis layer in order to stimulate and re-model collagen effectively by either RF microneedling mode(invasive) or RF matrix mode(non-invasive).

It can be used for both epidermal application and dermal applications in one platform without change any operational environment.



- Medical Grade Sterilization
- Gold Plate
- Single Use

### THE NEWEST GENERATION IN ADVANCE MICRONEEDLING

Micro-needle Fractional RF system is an ideal technology for skin lifting, tightening and rejuvenation by applying precisely controlled RF energy directly into various depths of skin with minimally invasive micro-needles.

It can be applied to all skin types, with very low risk of skin burns and PIH. This is a key difference when compared to optical based treatments(e.g. IPL/Laser).

F - 3\*3

- Eyes
- Nose Area

F - 5\*5

- Face Area

F - 9\*9


- Body Area



# INSULATED & NON INSULATED


## RF Target Area

### -Dermis

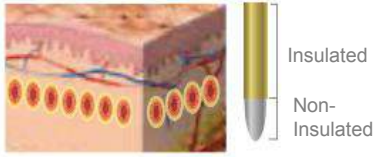

## RF Target Area

### -Epidermis & Dermis

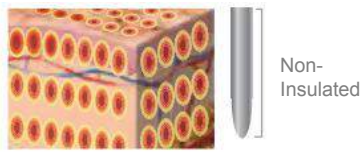

## Insulated


- Less painful
- Treat Dermis only
- Less redness, easy to recovery
- Less result

- Suitable for person who need go to work everyday

## Non-insulated

- Painful
- Treat Dermis and Epidermis
- More recovery time, easy scab
- Better result

- Suitable for person who want good result

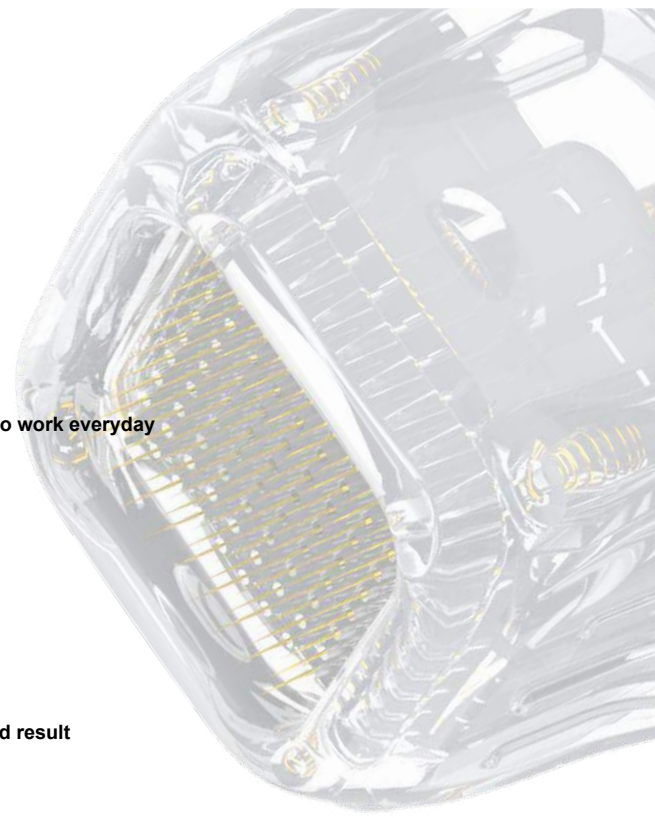

## APPLICATIONS



PORE



SCAR



LIFTING



WRINKLE





TIGHTENING



# ADVANTAGE-VACUUM

Pinxel-V update RF microneedling with vacuum, push the topical deeper with one kind dual air chamber design, air is captured through enhanced tissue engagement, microneedles creat channels and RF energy is delivered, as needles retract, air is released towards the skin, driving topicals deeper into the dermis, resulting in enhanced topical penetration up to 67% with the M and F Type needles.

Pinxel-V takes RF microneedling to new levels. The versatility this system provides allows you to customize treatments based on your patients' specific needs – all skin types, anywhere on the body, any time of the year.

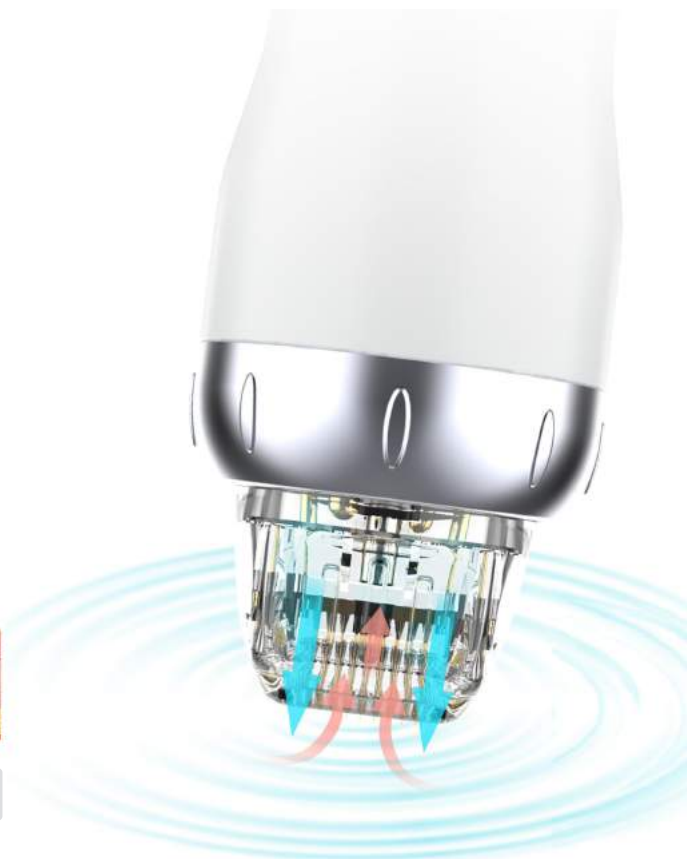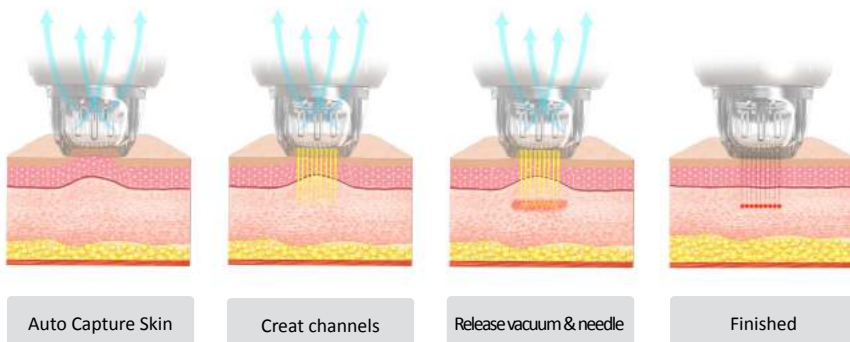




# MESO

MESO SKINJECTOR

## MESO SKINJECTOR


Water mesotherapy is to inject hyaluronic acid into skin. Hyaluronic acid contains moisture, 200-300 times as rich as one's body. So that it makes the skin smooth, soft, bright and hydrated after the injection. In order to keep the skin hydrated, ladies always choose varied skin care products. Actually, skin care products only hydrate moisture on the cutin. So that over smearing those products cannot deeply provide moisture for the skin. As we know, acne, wrinkle and inelasticity would appear if skin lacks moisture. Water mesotherapy will solves all the problems by reason of lacking moisture, such as sallow skin, color spot, wrinkle... And it will supply moisture for dermis. Therefore water mesotherapy will works well on whitening, moisturizing and renewing inelasticity.

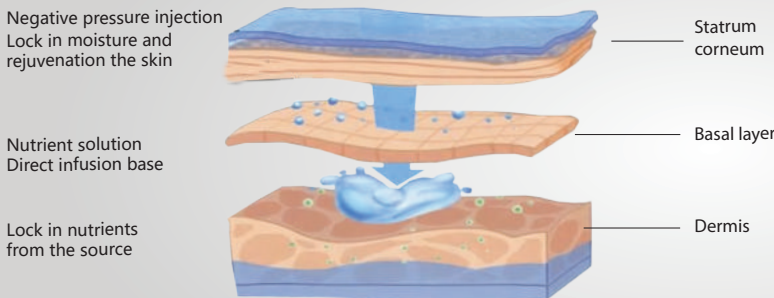


## TREATMENT PRINCIPLE

Micro-needle mesotherapy gun using voice vacuum technology, accurately provides nutrients needed for skin in the depth of 1.28mm under the dermis. It is the principle of lifting the skin slightly by using sound pressure before micro-needle entered. The needle injects nutrients (such as hyaluronic and collagen) by using multi-needle sticking in the accurate depth after keeping skin tight. Injection rate can be easily adjusted through the button on OLED panel behind the device. Further more, the pressure of syringe had been liberalized before micro-needle was out, it will not bring the loss of nutrients.

TECHNICAL SPECIFICATION	
Voltage	1 10V/240V 50Hz/60Hz
DA	12V/5A
Power	50W
Display	4.3 inches LCD touch screen
Multi-needle	3 needles, 4 needles, 5 needles, 9 needles
Mode	Low speed / High speed / 10-250 times/min
Advantage	Sing needle, Multi-needle, Vacuum, 3 in 1
Single needle caliber	0.25mm, Depth: 0.25~10mm
Multi-needle caliber	0.25mm, Depth: 0.25~5mm
Vacuum FloW	25L/min Pressure: 90Kpa

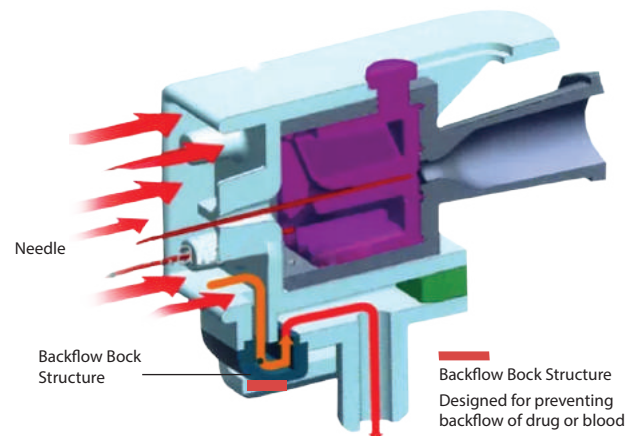




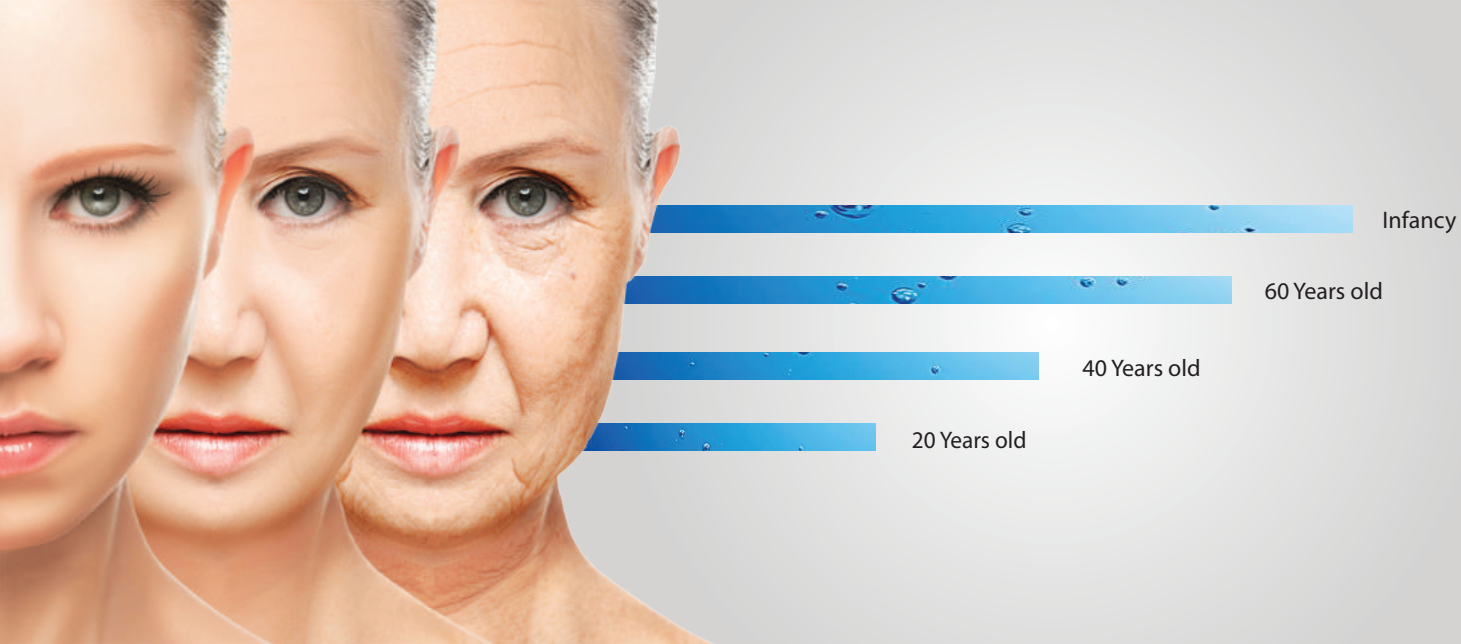


□ **ADVANTAGE**

- 1. Our handle is South Korea imported, as well as the needles.
- 2. It has 5 pins injection technology which can speed up the treatment process.
- 3. Adjustable injection depth, doses, speed and suction intensity.
- 4. Anti-back flow technology is applicative for medicine and blood.
- 5. Language: Chinese,English,Russian,French,German,Arabic.

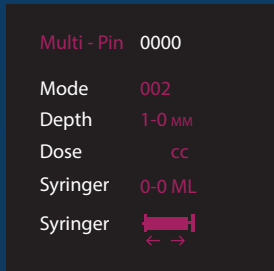

HANDLE DISPLAY

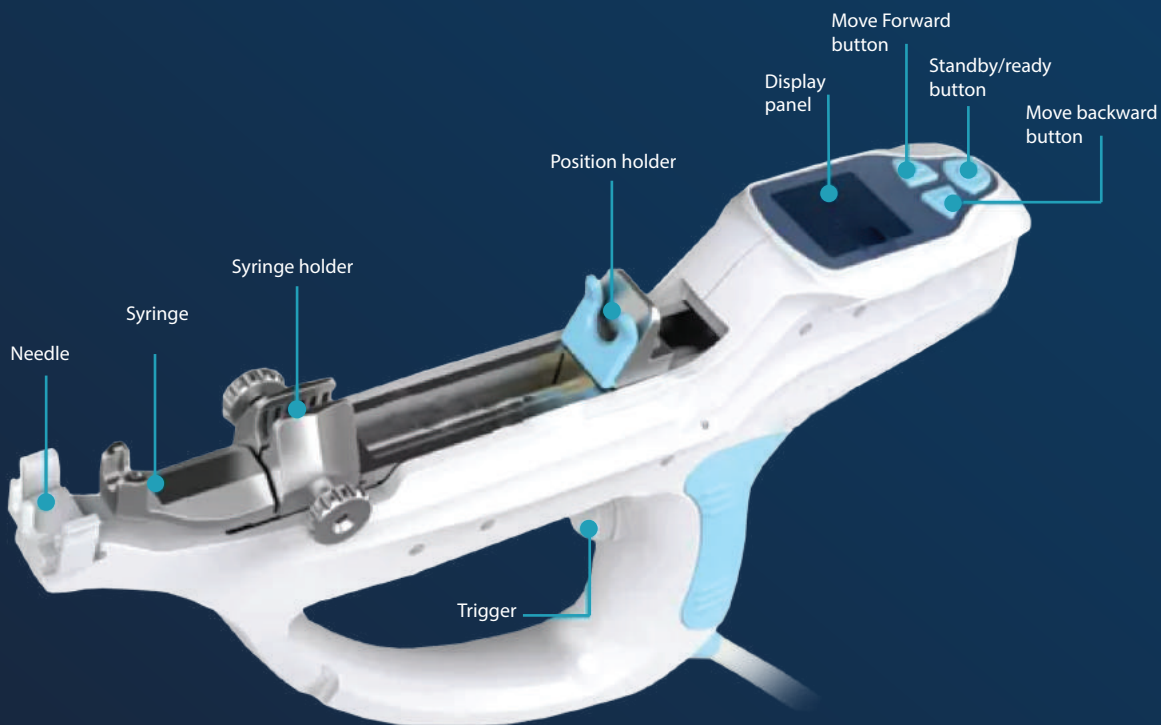

**THE EFFICACY OF NEEDLE MESOTHERAPY**

Micro-needle mesotherapy use of the minimally invasive ways to stimulate bone colleagen, the skin without damaging the epidermis.

- 1 Using the principle of negative pressure adsorption, the nutritional factors are instantly injected into the skin to replenish nutrients
- 2 The molecular particles are smaller than the pore diameter and can penetrate deep into the basal layer
- 3 Choose the right nutrient solution for different skin problems
- 4 Nutrients go directly under the skin, hydrating and locking in water



# FACECODE


Non-Invasive Skin Tightening Treatment



## HIFU FACECODE

Non-Invasive Skin Tightening  
Treatment

FACECODE A non-invasive treatment, goes deeper to stimulate production of your body's own collagen and elastin, for natural looking results without surgery or downtime.

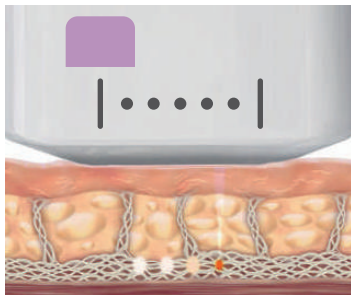

### HIFU treatment Mechanism

HIFU is delivered into SMAS(Superficial Muscular Aponeurosis System),Fascia, deep dermis and substaneous fibrous tissue and generates the temperatureincreased, when the temprature grows more than 47°C (the Threshold of protein denaturation), the coagulation zone is generated. the zone is filled out by the tissue/muscle surrounded and in this moment, the skin muscle gets elastic and the wrinkle is tightened.

#### Hand Piece

Thinner & Curved line.  
Very comfortable to handle  
Double switches.

#### Transducers

Super Thin

#### Display

LCD with touch screen

#### Main Body

Smiple & Smart Design



# THEORY


HIFU is a non-surgical radio frequency procedure that can rejuvenate your skin, help smooth, tighten and contour skin and decrease the sagging skin of aging on the face and body.

HIFU procedure stimulates collagen production by delivering energy to the skin's foundational layer typically addressed in cosmetic surgery-without cutting or disrupting the surface of the skin. Coagulation zone into 4.5mm Depth(SMS layer) and 3mm.

Make the higher temperature than 400°F, triggers a natural response under the skin, immediately tightening effect, jump starting the regenerative process that build collagen & elastin.



Mutifunctinal, a real low invest, high return machine.

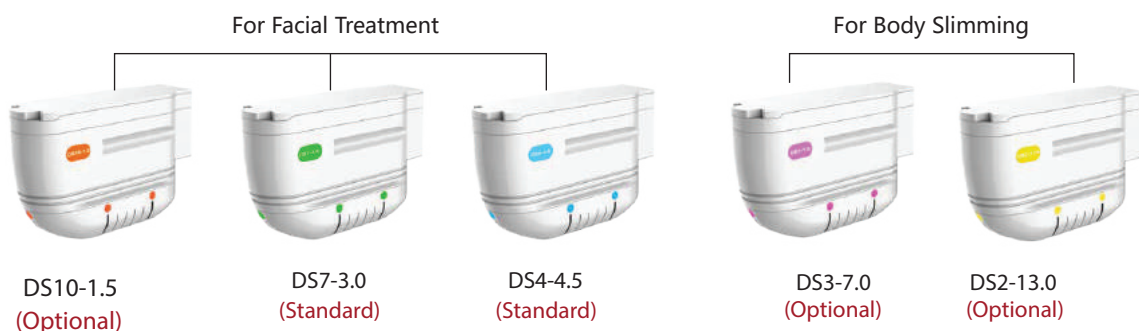

## Technical Specifications


Irradiation Type	High Intensity Focused Ultrasound
Frequency	4MHz, 7MHz
Max Energy	15 Inch True Color Touch LCD ADSS System
LCD Screen	2.0j/m2
Standard Cartridges	4MHz 4.5mm, 7MHz 3.0m---10,000 shots
Optional Cartridges	2MHz 13mm(Body/Reduces Fat)
	3MHz 7mm(Body/Reduces Fat)
	10MHz 1.5mm(Eye Area)
Length	Face 5~25mm(1mm step)
	Body 5~25mm(7.0mm, 13mm)
Voltage	AC110-220V 50/60Hz

## APPLICATION

### Suitable for what kind of Conditions ?

Treat sagging skin, aging phenomenon for any areas such as forehead, eye brow, cheekbones, apples muscle, cheeks, neck and body

- Skin tightening and rejuvenation
- Non-Invasive Lifting for Face
- Body Contouring and tightening

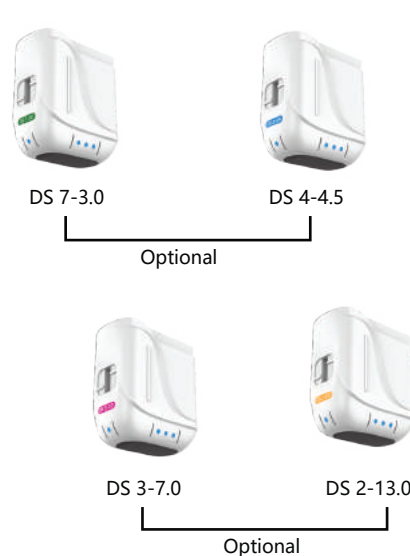



# Visage

The Latest Facial Treatment To Revitalize  
Skin in a Non-Invasive Way

**Visage** targets & lifts frontalis muscles, corners of the mouth, zygomaticus major and minor muscles, matching the contours of every patient's face, activating & contracting these individual muscles. Meanwhile, special temperatures induce new collagen & elastin growth to resist aging. The combined effects restore and elevate the support of facial tissues, leading to an overall improvement in facial lifting.

EM Time

The first & only facial procedure to treat:

Needle-Free  
Solution

Simultaneous Skin  
&  
Muscle Enhancement

Natural  
Lift

Effective Wrinkle  
Reduction

Full Face  
In 20 Minutes

Reveal your beauty, naturally.

# Simultaneous Skin & Muscle Enhancement



Injection-Free Solution



Hands-Free



Quick Results



No Downtime



- One course offers long-term results with the overall improvement of the structure of the face and skin quality
- No downtime required, nor the need for at-home recovery
- Gentle and pain-free non-invasive facial treatment requires no scalpels, no needles, no fillers
- It is safe for most people, including those with sensitive skin
- Equipped with unique technologies to target all facial layers, delivering non-invasive facial rejuvenation

Reveal your beauty, naturally.

## 2 Technologies In 1 Treatment

### Muscles & Skin

Visage works on the framework of the face from the inside out and targets facial layers with complicated cutting-edge tech. It restores and elevates the support of facial tissues by selectively stimulating muscles and increasing the density and quality of muscle structure, enhancing the facial structure, and simultaneously remodels and smooths the skin by increasing collagen and elastin fibers, which results in a reduction of wrinkles and fine lines while tightening the skin. With EM Time treatment, you can enjoy long-lasting benefits and a rejuvenated, lifted, toned, and strengthened appearance.



# WAVELENGTH PRO-X



Safety



Painless



Simple  
Operation



Hair Removal

## INTRODUCTION

Diode laser hair removal machine is a laser that produces a larger beam of highly concentrated light. The laser beam is well absorbed by the pigment located in the hair follicles. During the procedure, the laser pulses for a fraction of a second, allowing the hair to absorb the energy and heat up. As it heats up, the hair shaft and bulb are damaged which significantly impedes its ability to re-grow. And with the large spot size, the process treats numerous hair follicles simultaneously, making the treatment much faster than traditional IPL and other hair removal alternatives.

## Technical specifications

### Basic Operation

Laser type	Diode Laser
Laser Bars Power	1600W (160w x 10pcs)
Wavelength	755nm+808nm+1064nm

### Technical characteristics

Display	15" Color Touch LCD Display
Spot Size	12x 20 mm
Pulse width	2-150ms adjustable
Output power	300W
Energy	2-100 j / cm2 adjustable
Frequency	1-10 Hz
Voltage	AC 220V/50Hz
Cooling mode	Close- Cycle water cooling
Cooling	Water+air+semiconductor
Temperature of Probe	0~30 C




## 3 Combined Wavelengths Larger Spot Size

755nm	Suitable for People with fair skin, good effect on eyebrow and upper lip.
808nm	Suitable For Yellow skin or light skin.
1064nm	Suitable for dark skin.

## THEORY

The diode laser system can emit laser light of 808nm wavelength. It can be absorbed by the melanin located in the hair follicle so as to destroy the hair follicle. The hair shaft can also absorb the energy from laser light, so as to be destroyed and gasified. then the unwanted hair can be removed completely and permanently.

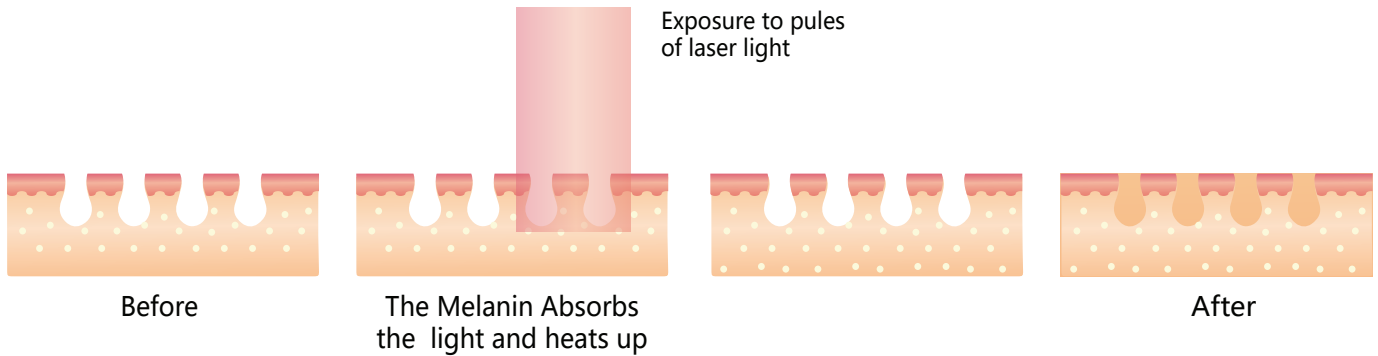

## 3 Three in one 3wave

### 3 Wave Technology

Featuring its exclusive 3 wave technology, wavelength pro-x combines the three most effective laser wavelengths into a single applicator, simultaneously targeting different tissue depths and anatomical structures within the hair follicle. By combining the absorption and penetration levels of different wavelengths, along with extended treatment coverage, increased comfort and low maintenance requirements, wavelength pro-x is the safest and most efficient hair removal treatment available today 3 Dimension technology for one superior hair removal solution.



#### 755 nm Alexandrite Laser

Powerful absorption, combined with continuous cooling, Treatment coverage, comfort, and low maintenance, Enables the treatment of a wide range of hair types and colors, Especially light-colored and thin hair.



#### 808 nm Diode Laser

Complete more treatment sessions in less time with this large spot size applicator. Treat larger areas faster, while offering patients a safe and highly effective hair removal procedure.



#### 1064 nm Diode-Yag Laser

The YAG 1064 uses a diode laser based on an ND : YAG wavelength, allowing for moderate energy absorption by the melanin chromosphere, and is optimally suited for patients with darker skin tones.






# CODE

A é s t h é t i c s

CERTIFIED BY

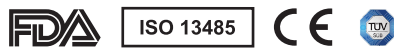

---

[www.codeesthetics.com](http://www.codeesthetics.com)

---

## CODE AESTHETICS

CORPORATE OFFICE  
B2-34/A First Floor  
Yamuna Vihar Delhi-110053

CONTACT :  
+91 9716830460  
Info@codeesthetics.com

